

Mixing Cells Protocol

CELLINK A Series

This is a suggested procedure, please adjust according to your experimental needs.

Protocol aim

The aim of this protocol is to provide instructions for how to mix cells with any bioink within the CELLINK® A Series, both small volumes below 2 ml and large volumes up to 10 ml. The CELLINK® A Series includes CELLINK® A and CELLINK®. Both bioinks can be purchased with the STARTINK-Kit.

Material needed

- Cells in suspension
- Cell culture medium
- CELLINK® A or CELLINK® A-RGD*
- Pipette and pipette tips
- 3 ml syringes with luer lock connections
- Female/female luer lock adaptor*
or
- CELLMIXER*

1

* The product can be purchased in the CELLINK shop at www.cellink.com/store/.

Protocol

This protocol is adjusted for mixing either 1 ml or 3 ml of bioink with cells to a final cell concentration of 10 million cells/ml bioink. For other quantities and cell concentrations recalculations need to be made but the same protocol can be followed.

Step	Title	Material	Description
1	Prepare cell suspension	<ul style="list-style-type: none">- Cells in suspension- Culture medium	<p>If preparing for quantities < 2 ml of CELLINK® A bioink.</p> <ul style="list-style-type: none">- Resuspend 11 million cells in 100 µl cell culture medium if mixing with 1 ml bioink.- Move on to Step 2a. <p>If preparing for quantities > 2 ml of CELLINK® A bioink.</p>

			<ul style="list-style-type: none"> - Resuspend 33 million cells in 300 µl cell culture medium if mixing with 3 ml bioink. - Move on to Step 2b.
2a	Mixing small volumes CELLINK® A variation with cells	<ul style="list-style-type: none"> - 1 ml CELLINK® A variation - Cells in suspension - 3 ml syringes - Female/female luer adapters 	<p>At this point, mix ten parts bioink with one part cell suspension, taking care not to introduce air bubbles to the mixture.</p> <ul style="list-style-type: none"> - Transfer the 100 µl cell suspension to a 3 ml syringe using a female/female luer lock adaptor. - Transfer 1 ml of bioink to a 3 ml syringe using a female/female luer lock adaptor. - Attach the bioink syringe to the syringe with cell suspension. - Carefully mix the bioink with the cell suspension by gently pushing the bioink back and forth between the syringes. - Transfer the cell containing bioink back to the cartridge and cap it. - <i>Video link for a detailed illustration on how to perform the mixing process:</i> https://www.youtube.com/watch?v=NmdOTNLrV-Q <p>Note: To avoid an air gap when mixing the bioink and the cell suspension, carefully pre-fill the luer lock adaptor with CELLINK® A bioink before attaching the syringe with the cell suspension.</p>
2b	Mixing larger volumes CELLINK® A with cells	<ul style="list-style-type: none"> - 3 ml CELLINK® A variation - Cells in suspension - CELLMIXER 	<p>At this point, mix ten parts bioink with one part cell suspension, taking care not to introduce air bubbles to the mixture.</p> <ul style="list-style-type: none"> - Transfer the 300 µl cell suspension to the 1 ml cell syringe (PART 1) using a female/female luer lock adaptor. - Transfer 3 ml of bioink to the 12 ml syringe (PART 2) using a female/female luer lock adaptor. - Clip both syringes to the Dispensing unit (PART 3). - Connect the two syringes to the Mixing unit (PART 4), then connect the Empty cartridge (PART 5) to the Mixing units other side. - Apply gentle pressure onto the Dispensing unit to mix the content of both syringes into the empty cartridge.

			<p>- Video link for a detailed illustration on how to perform the mixing process using the CELLMIXER: https://www.youtube.com/watch?v=CmSYL1-olt/</p> <p>Note: PART 2 is prefilled with 3 ml of bioink if using the STARTINK-Kit.</p>
--	--	--	---

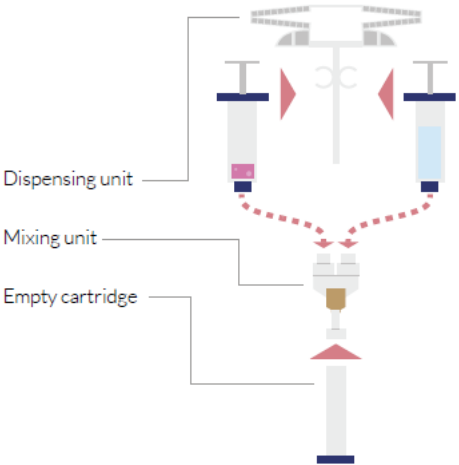


Figure 1. Illustration of how to assemble the CELLMIXER components.